
waterpik®

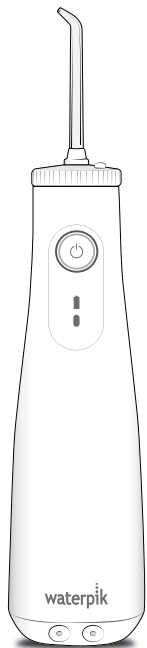Waterpik® Water Flosser Model WF-10 Series

www.waterpik.com

IMPORTANT SAFEGUARDS

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING.

DANGER:

To reduce the risk of electrocution:

- Always unplug product after using.
- Do not handle plug with wet hands.
- Do not place in or drop into water or other liquid.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub, or sink.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Check the power supply cord for damage before the first use and during the life of the product.

WARNING:

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not charge this unit with any charging devices other than the one provided.
- If product is opened/disassembled for ANY reason, warranty is **VOID**.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into any water. Contact Waterpik International, Inc. or your local distributor which can be found at www.waterpik.com/intl.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Only use tips and accessories recommended by Water Pik, Inc.
- Do not drop or insert any foreign object into any opening or hose.

- Keep electrical cord away from heated surfaces.
 - Magnetic field: Pacemakers and similar implanted medical devices may be affected by magnet.
 - Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
 - Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. Use of these can reduce product performance and will shorten the life of the product. Use of these products may void warranty.
 - Remove any oral jewelry prior to use of this product.
 - Do not use if you have an open wound on your tongue or in your mouth.
 - If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
 - If the Pik Pocket™ tip (not included in all models) separates from the shaft for any reason, discard tip and shaft and replace with new Pik Pocket™ tip.
 - Do not use this product for more than 5 minutes in each two-hour period.
- Closely instruct and supervise children 8 years and above and individuals with special needs in the proper use of this product.
 - Children should be supervised to ensure that they do not play with the appliance.
 - There are no consumer serviceable parts in this appliance and it does not require maintenance.
 - Fill reservoir with warm water or other dental professional recommended solution only.

SAVE THESE INSTRUCTIONS.

TABLE OF CONTENTS


Product Description	5
Getting Started	6
Cleaning and Troubleshooting Your Waterpik® Water Flosser	9
Limited Two-Year Warranty	10

1. Tip rotator
2. Power button
3. Mode/pressure button
4. Battery charge indicator
5. Tip
6. Tip eject button
8. Reservoir fill door
9. Reservoir

Charger


- 20024431-3 White
- 20024431-4 Black

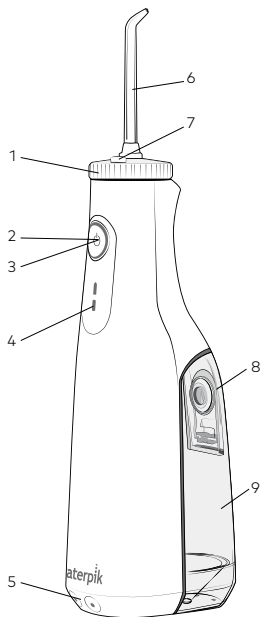TIPS

- | | |
|---|--|
| <p>A. Orthodontic Tip*
- Braces/General Use</p> <p>B. Pik Pocket™ Tip*
- Periodontal Pockets/Furcations</p> <p>C. Plaque Seeker™ Tip*
- Implants/Crowns/Bridges/
Retainers/General Use</p> | <p>D. Toothbrush Tip*
- General Use</p> <p>E. Classic Jet Tip**
- General Use</p> <p>F. Tongue Cleaner**
- Fresher Breath</p> |
|---|--|





NOT ALL TIPS INCLUDED IN ALL MODELS

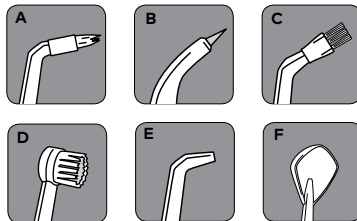*Replace Every 3 months **Replace Every 6 months

Replacement tips can be purchased from your local distributor, which can be found at www.waterpik.com/intl.

GETTING STARTED

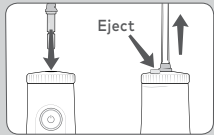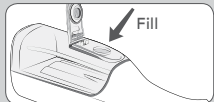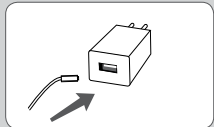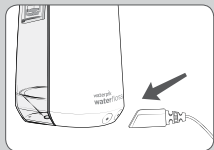Charging Unit Prior to First Use

Plug the power cord into the wall and place the magnetic charger on the front of the unit for approximately 4 hours (see image). Ensure contacts are clean/dry before connecting the charger. The LED light will display a single, rapid blink every 3 seconds when unit is fully charged and the charger is connected. A full charge will last for approximately 1 week of typical use. When the batteries are in need of charging (approximately 30% charge left) the charge indicator light will blink 8 times quickly after use to indicate that you should recharge the unit. A complete recharge will take 3-4 hours. If you use your Water Flosser once a day or less, it is recommended you charge it once a week. If you use it twice a day or more, charge the unit more frequently.


Filling the Reservoir

Remove the magnetic charger from the unit prior to filling reservoir and product use. Lift flip top on reservoir and fill the reservoir with lukewarm water.

Inserting and Removing Tips

Insert tip into the center of the knob at the top of the Water Flosser handle. Press firmly until the tip clicks into place. To remove the tip from the handle, press the tip eject button and pull tip from the handle.

CAUTION: Do not press tip eject button while unit is running.

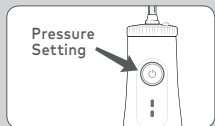Adjusting the Pressure Setting

- First press (from OFF) starts unit on "LOW".
- Second press turns unit to "HIGH".
- Third press turns unit "OFF".

The single LED exhibits the following behavior:

- Solid - If motor is running and battery is greater than half full.
- Single blink every 2 seconds - When charger is connected and battery is fully charged.

Recommended Technique

Lean low over sink and place tip in mouth. Aim the tip toward gumline, and turn the unit ON (I).

For best results, start with the back teeth and work toward the front teeth. Glide tip along gumline and pause briefly between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth.


Direct the jet stream at a 90-degree angle to your gumline. **Slightly close lips to avoid splashing** but allow water to flow freely from mouth into the sink. Keep unit upright during use for best results.

When Finished

Turn the unit OFF. Empty any liquid left in the reservoir.

If desired, the reservoir may be removed from the power handle by simply pulling down on the reservoir to remove.

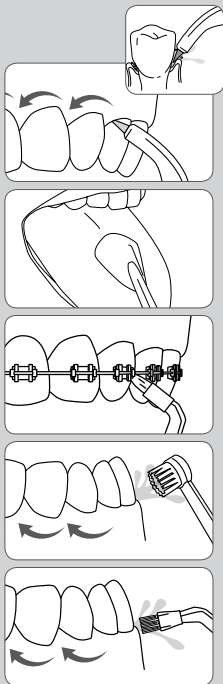NOTE: NOT ALL TIPS INCLUDED WITH ALL MODELS.

TIP USAGE

Pik Pocket™ Tip

The Pik Pocket™ tip is specifically designed to deliver water or anti-bacterial solutions deep into periodontal pockets.

To use the Pik Pocket™ tip, **set the unit to the lowest pressure setting**. Place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket.

Tongue Cleaner Tip

Set to lowest pressure setting, place tip in the center/middle of your tongue about half way back. Pull forward with light pressure. Increase pressure as you prefer.

Orthodontic Tip

To use the orthodontic tip, gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth.

Toothbrush Tip

Place the toothbrush tip in mouth with the brush head on the tooth at the gumline. The toothbrush tip can be used with or without toothpaste. Using a **light pressure** (bristles should not bend) massage the brush back and forth with very short strokes – much like you would with a manual toothbrush.

Plaque Seeker™ Tip

To use the Plaque Seeker™ tip, place the tip close to the teeth so that the bristles are gently touching the teeth. Gently glide tip along gumline, pausing briefly between teeth to gently brush and allow the water to flow between the teeth.

Using Mouthwash and Other Solutions

Your Waterpik® Water Flosser can be used to deliver mouthwash and antibacterial solutions. For ease of adding other solutions remove the reservoir (see "When Finished" section) and add from top opening on reservoir. If you are adding additional water you may replace the reservoir and fill reservoir fully with the fill door. After using any special solution, rinse unit to prevent clogging by partially filling reservoir with warm water and running unit with tip pointed into sink until unit is empty.

Cleaning

Never immerse unit in water. Before cleaning, unplug from the electrical outlet. Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser to wipe exterior.

The water reservoir is removable for easy cleaning and is top rack dishwasher safe.

It is recommended that the reservoir is removed and cleaned weekly.

Removing Hard Water Deposits / General Cleaning

Hard water deposits may build up in your unit, depending on the mineral content of your water. If left unattended, it can hinder performance. Cleaning internal parts: add 1 tablespoon of white vinegar to a full reservoir of warm water. Point the handle and tip into sink. Turn unit ON and run until reservoir is empty. Rinse by repeating with a full reservoir of clean warm water.

This process should be done every 1 to 3 months to ensure optimal performance.

Battery Care and Disposal

If the Water Flosser is to be stored for an extended period of time (six months or more), be sure to fully charge it prior to storage. The unit contains a Nickel Metal Hydride battery. At the end of the battery life, please recycle the unit at your nearest recycle center per local or state requirements. The battery is not replaceable.

LIMITED TWO-YEAR WARRANTY

Service Maintenance

Waterpik® Water Flossers have no consumer-serviceable electrical items and do not require routine service maintenance. For all your warranty and accessory needs, please email InternationalService@Waterpik.com. Please do not ship product to Waterpik International, Inc. offices, this will delay service.

Refer to the serial and model numbers in all correspondence. These numbers are located on the back of the handle and bottom of the reservoir respectively.

Still have questions? Go to www.waterpik.com/intl to find your local distributor.

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 2 years from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will replace at its discretion the product, which in its opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. Repair or replacement services provided for by the guarantee will be carried out by the local distributor, whose details may be found on www.waterpik.com/intl. Upon request of the consumer, a copy of the guarantee complete with the local distributor's information will be provided by email or mail. This limited warranty excludes accessories or consumable parts such as tips, etc.

This warranty gives you specific legal rights. You may also have other rights that vary depending on the law in your place of residence. To locate your local distributor, go to www.waterpik.com/intl.

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- Do not recharge the battery of this product with any other charging device than the one provided.
- If the product is opened or disassembled for ANY reason, the warranty is **ANNULÉE**.
- Do not use this device if the plug or the power cord is damaged, if it does not function properly, if it has fallen or been crushed, or if it has fallen into water. Please contact Water Pik, Inc. or your local distributor at www.waterpik.com/intl.
- Do not direct the water jet of the hydrojet under the tongue, in the ear, in the nose or towards any other delicate area of the body. This device is capable of producing pressures that could cause serious injury in these areas.
- Use the device in the manner indicated in the user manual or follow the recommendations of your dentist.
- Use only the tips, toothbrushes and accessories recommended by Water Pik, Inc.

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©2019 Water Pik, Inc.
Church & Dwight Co., Inc.
Waterpik, International, Inc.
Dwight Co., Inc.